

Medipoint® 125 Central Alarm

The BeaconMedæS Medipoint® 125 Medical Gas Central Alarm system is designed to provide full monitoring of the complete medical gas system installation. An individual Medipoint 125 alarm panel displays up to 5 gas sources in the 'Normal' and 4 'stage' alarm conditions, or may be used to display up to 20 individual point alarms. Line contact monitoring circuits constantly verify the integrity of input sensors and interconnecting wiring circuits.

Antimicrobial Fascia

The fascia is made from Autotex AM. Microban® technology is incorporated into the Autotex textured hardcoat during the manufacturing process. When bacteria come into contact with Autotex AM with Microban®, the antimicrobial function disrupts the bacterial cell wall, killing or inhibiting bacterial growth.



Autotex AM has been tested to AATCC Test Method 100⁷, and exhibited biocidal passes against a wide range of microbes including MRSA and E. coli.

Surface Mount Technology

The Medipoint 125 now comes with advanced printed circuit boards using surface mount technology as standard. This allows cooler operating temperatures to be maintained, increasing the reliability and substantially lengthening service life.

*Microban and the Microban symbol are trademarks of the Microban Products Company, Huntersville, NC, USA



*Life
is in the
details.*



Expands with your Hospital


Features like the pre-punched 1st fix box, interchangeable gas service legends and simple switch programming facility allow for system expansions to be accommodated without the need to replace or add anything.

Simple Reprogramming

The Medipoint® 125 Alarm is easy to reprogramme on site if modifications are made to your medical gas system. Simple DIL and rotary switches are used to set the panel ID, column inputs and the mute facility - no mess, no fuss - simple.

Adjustable Volume Control

The alarm also includes a volume control so that the volume can be set to suit the area of installation. This feature is of great benefit in areas where alarms are positioned close to monitored nurse bases where patients should not be disturbed.

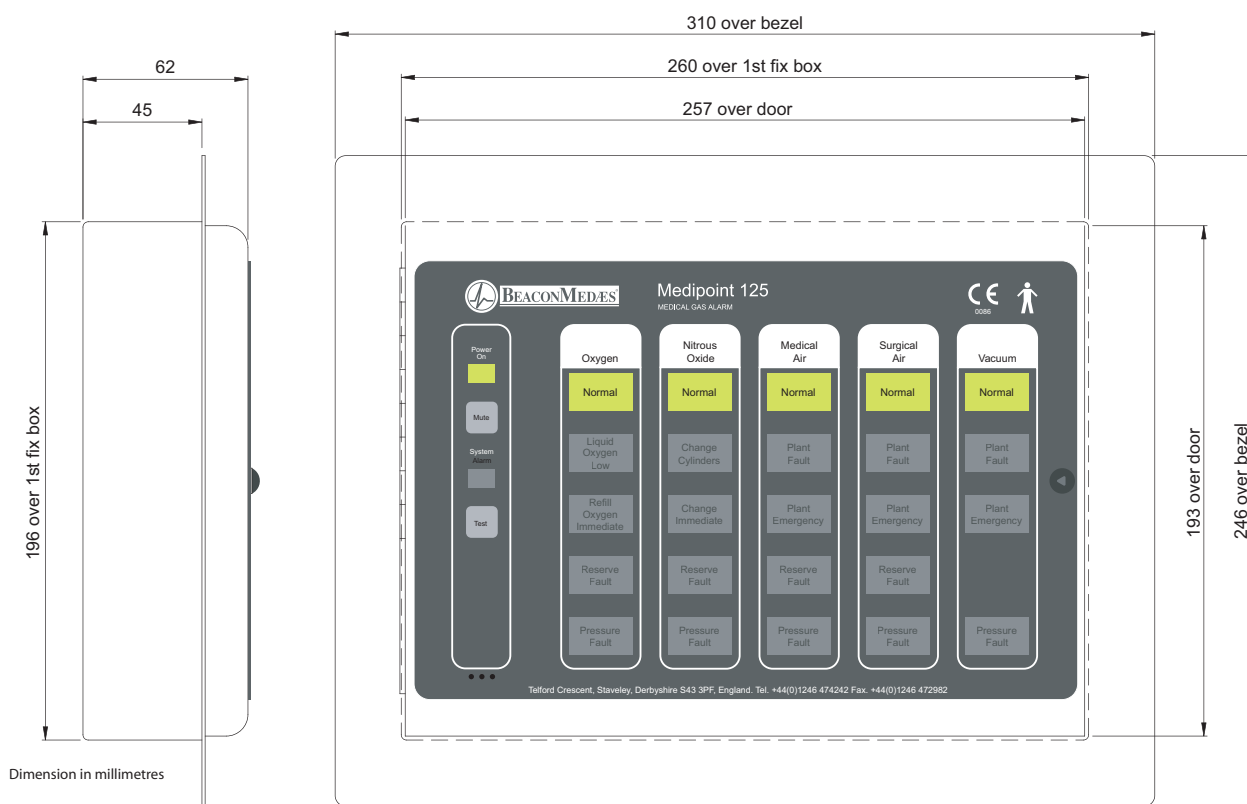

Medipoint® 125 Central Alarm

Description	Part Number	Description	Part Number
1st fix box	1826561	2 pair data cable (<500m)	1824637
Bezel for flush installations	1826564	2 pair data cable (>500 <1000m)	1829807
2nd fix for 1-5 columns - 230V	1829815	6 core screened input cable	1829806
2nd fix for 1-5 columns - 110V	1829820	Remote Audible	1826610
Relay interface 2nd fix - 230V	1829808	Signal booster (IP67)	2004133
Relay interface 2nd fix - 110V	1829809		

BeaconMedæS
A company within the Atlas Copco Group
Telford Crescent
Staveley, Derbyshire
S43 3PF, England

Tel: +44 (0) 1246 474 242
Fax: +44 (0) 1246 472 982
Email: gbn.info@beaconmedaes.com
Web: www.medaes.com

Medipoint® 125 Central Alarm

